



High Flow Heated Respiratory

Humidifier


This product is suitable for patients with spontaneous breathing, realizing treatment by providing heated and humidified gas of a certain flow. These patients include humidification treatment, oxygen therapy, endotracheal intubation and tracheostomy. This product cannot be used for life support.

Flow rate

The stable flow rate can be accurately adjusted, according to the patient's individual ventilation demand;

Low Flow Mode 2-25 L/Min (pediatric mode);

High Flow Mode 10-80 L/Min (adult mode);

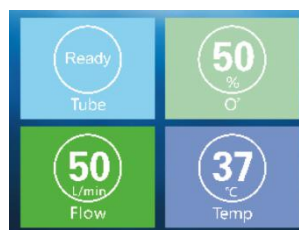

Monitoring

The display provides information on Current flow rate,

Temperature,

Therapy air FiO₂,

Mode of operation.

Oxygen


O₂ addition can be titrated /auto-control according to the patient's O₂ deficit.

Supply from any external O₂ source Up to 80 L/min.

Humidification

The level of humidification and the temperature of therapy air can be adjusted by the patient according to comfort.

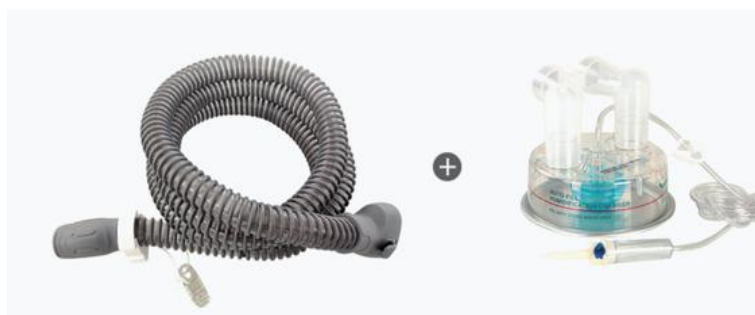

Dew point 31°C, 34°C, 37°C






PROMED HEALTHCARE LIMITED

Add: Shenzhen | Hongkong | China

Email: INFO@LLKMEDICAL.COM

| Item | HUMID-BM | HUMID-BH |
|--------------------------|--|--|
| PIC |  |  |
| Power supply | 220VAC+ 22V,50+1Hz | |
| Flow setting range | Low Flow Mode : 2-25 L/Min temp: 34°C High Flow Mode : 10-80 L/Min | |
| Temperature Control | 31°C ,34°C,37°C | |
| Oxygen % Measure | 21%-100% | |
| FiO2 Control | Manual titration control | Auto proportional valve control |
| Alarm Parameters | Low oxygen concentration alarm High oxygen concentration alarm Ambient temp alarm Choke alarm Alarm recovery | |
| Recommended Ambient Temp | 18°C ~28°C | |

| Model | Nasal cannula for adult | | | Nasal cannula for pediatrics | | Tracheostomy adapter |
|--------------|---|-------------|-----------|--|----------------|---|
| |  | | |  | |  |
| size | s | m | l | Infant size | Pediatric size | For Adult (flow rate :10-60 Lpm) |
| Cannula size | 3x1.5 mm OD | 6X4.8 mm OD | 5.4 mm OD | 2.9 mm OD | 3.2 mm OD | |
| Flow rate | 50 Lpm | 70 Lpm | 80Lpm | 20 Lpm Max | 25 Lpm Max | |
| package | 10 Pcs /box : 230x130x160 mm 70 Pcs /box : 530x230x290 mm | | | Sticker material :Hydrocolloid | | |

PROMED HEALTHCARE LIMITED

Add: Shenzhen | Hongkong | China

Email: INFO@LLKMEDICAL.COM