


Clinical Automatic BP Monitor DBP-01P

Features

- Pulse wave anti-interference technology, reading Accurate;
- Measurement during intelligent inflation, fast & comfortable
- Big LED display, convenient for user readiness.
- Indication of air pressure and temperature.
- Automatic reminding of elbow position upon elbow position sensor.
- Cuff cylinder is adjustable at up-down right-left angle, easy fits user body position.
- Date & time settable.
- Various memory & connection port: USB, RS232, ethernet, SD.
- Automatic BP level reminding on each reading upon hypertension standard of WHO.
- Modular design, easy maintenance.

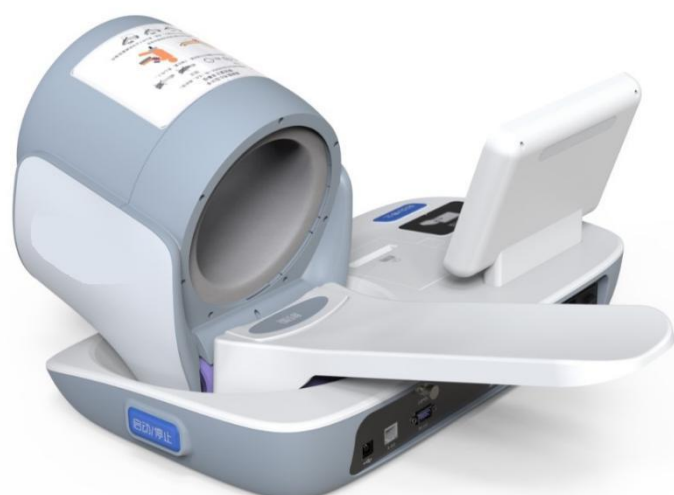

Specification

Technology	Intelligent Inflation Oscillometric
Display	8inch LED
Measurement Range	Systolic:40-260 mmHg Diastolic: 20-210 mmHg Pulse Rate: 40-200 bpm
Accuracy	±3 mmHg
Calibration	At least once per two years
Safety	Valve open automatically for deflation upon max inflation pressure 290mmHg. Pushing "Stop" button upon emergency, deflating quickly.
Memory	9999 readings
Power	110-240V AC, 50/60Hz

PROMED HEALTHCARE LIMITED

Add: Shenzhen | Hongkong | China

Contact: Dominic Qiu; Email: SALES55@LLKMEDICAL.COM_

PROMED 普尔迈

Printer	Thermal printing, Paper size 58mm
Cuff Circumference	17cm-42cm
Voice	Voice introduction of operation & reading, Language Option: Chinese, English. On/Off Option
Electric Shock Protection	Class I
Working Condition	Tem:+5°C -+40°C, RH: 10%-95%, Air Pressure 70-106Kpa.
Store Condition	Tem: -2°C -+55°C, RH: ≤95%, Air Pressure 70-106Kpa



Packing Detail

Unit size: 310mm x 476mm x 300mm

Weight: 6.5kg

PROMED HEALTHCARE LIMITED

Add: Shenzhen | Hongkong | China

Contact: Dominic Qiu; **Email:** SALES55@LLKMEDICAL.COM_

PROMED 普尔迈